

3M™ True Definition Scanner


Intraoral Workflow Guide



GLIDEWELL LABORATORIES

800-887-3580 • www.glidewelldental.com


How to Fill Out an Rx

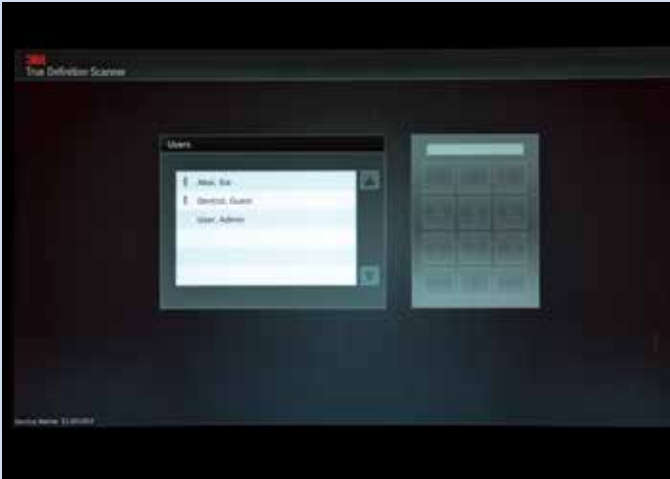

Indicate user, then provide user ID code



Press "Patients"

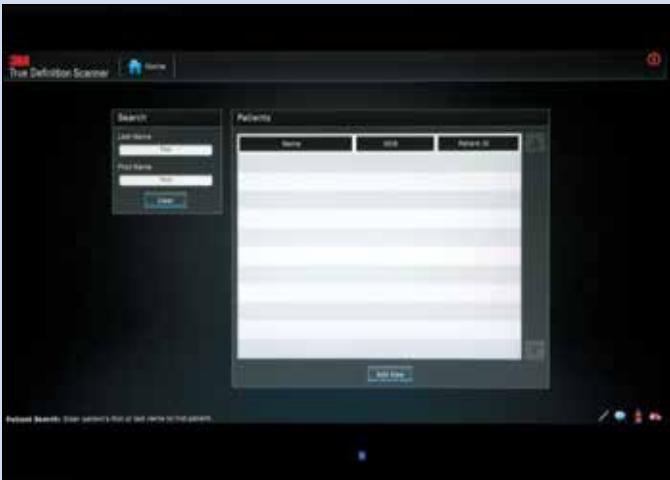

Enter the patient name, then press "Add New"




Fill in "Date of Birth" and any additional information, then press "Scan Now"

How to Prep Scan




Press appropriate arch and begin scanning

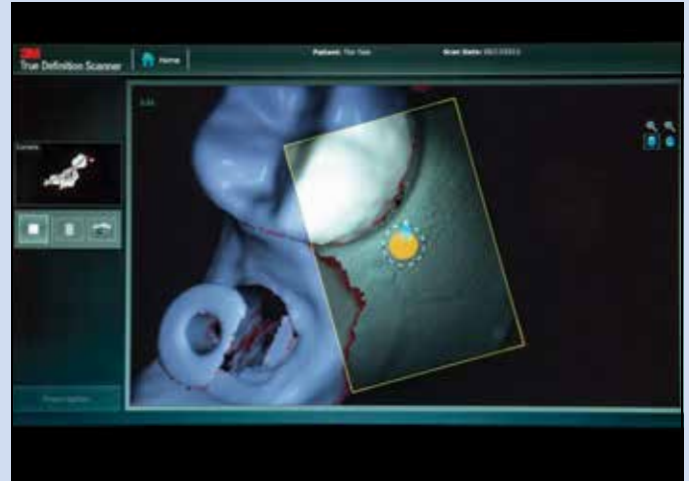

Scanning



When scan is complete, press the "add prep icon"

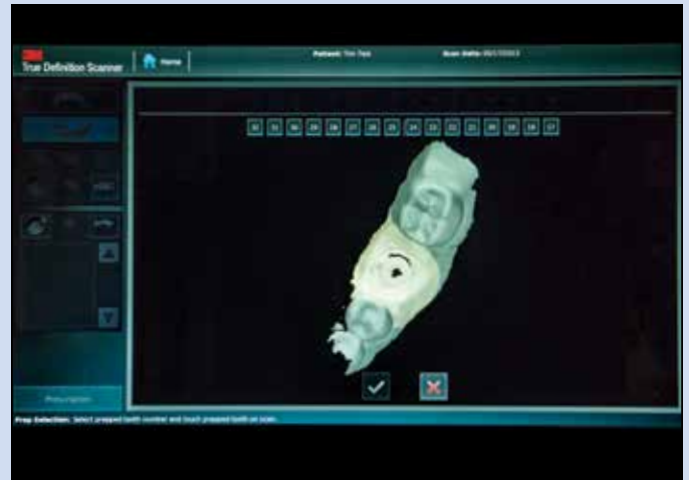


1. Touch and highlight the center of the prep. Prep and contacts will change color indicating high resolution
2. Press tooth number
3. Press green "✓"

How to Prep Scan



Review and accept by pressing "✓"

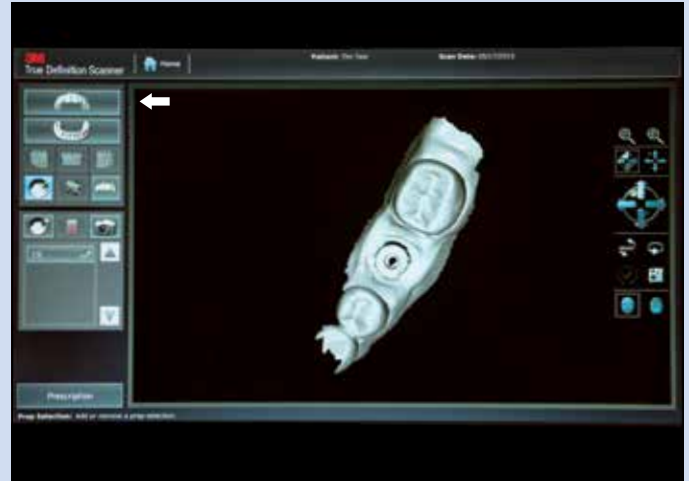

Press opposing arch and begin scanning



When scan is complete press "✓"




Press the appropriate bite icon

How to Prep Scan



Scan the bite

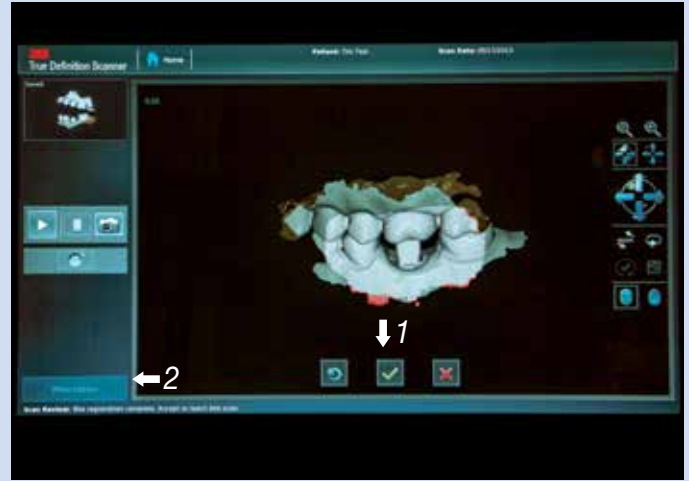

1. When bite is aligned, press the green "✓"
2. Press "Prescription"

How to Send Your 3M True Definition Scanner Files to Glidewell Laboratories




Select date request and press "Ok"



Select Glidewell Laboratories



1. Fill in notes with:
 - Implant system
 - Implant diameter
 - Restoration
2. Press "Touch to Select Restoration"

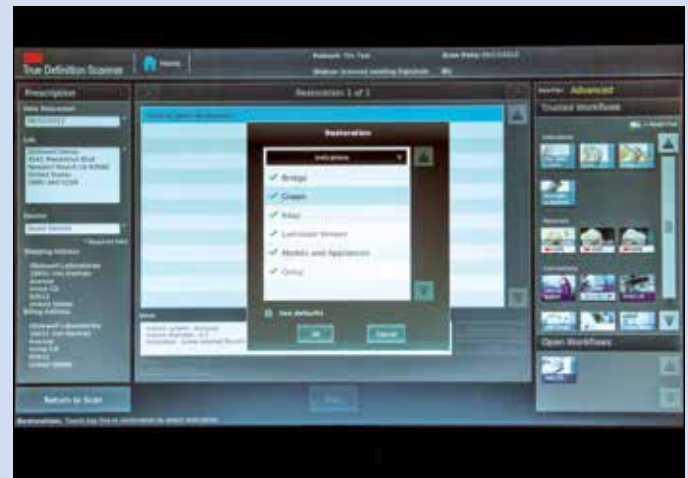

Press "Crown" then press "Ok"

How to Send Your 3M True Definition Scanner Files to Glidewell Laboratories



Highlight the restoration, then press "Ok"



Select the material from the right and press "Sign"



Enter the PIN, then press "Sign"



Press "Ok"



GLIDEWELL LABORATORIES

800-887-3580 • www.glidewelldental.com

GL-3390-0613