

3Shape TRIOS[®]

Intraoral Workflow Guide



GLIDEWELL LABORATORIES

800-887-3580 • www.glidewelldental.com

How to Fill Out an Rx




Select appointment time



1. Select Glidewell Lab
2. Fill out the patient information, fill out implant information
3. Select the scanner icon

How to Prep Scan

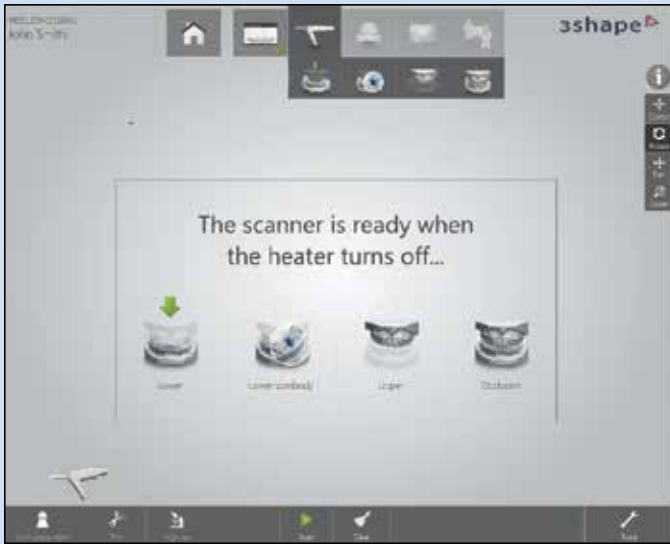

1. Select "Prep Arch"
2. Begin scanning



When scan is complete, press tooth #19 icon then press the scan abutment



Scan abutment indicated




1. Select scan body icon
2. Void will appear on the model
3. Once again, scan void area – this will be in high resolution

How to Prep Scan




When scan is complete, press #19 icon then press the scan abutment

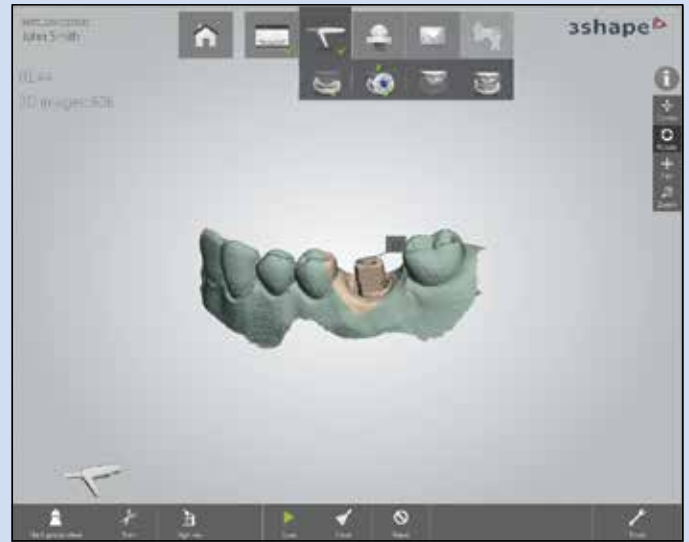

Scan abutment indicated



Select "Upper"

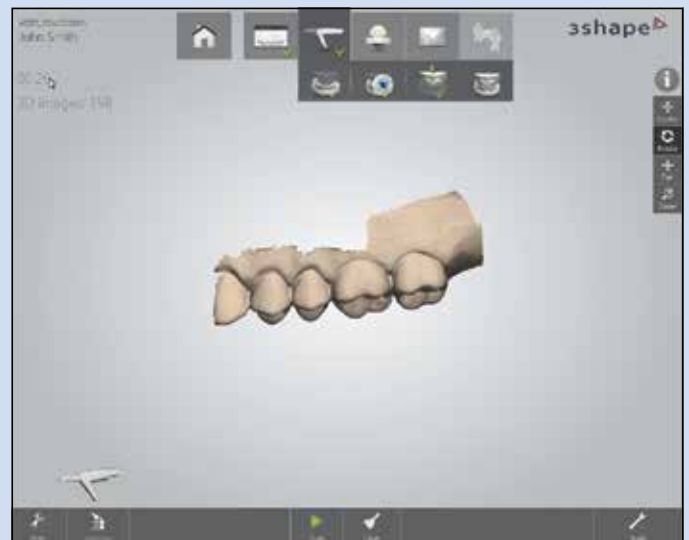


Flip the scan tip over 180° then scan the opposing arch


How to Prep Scan



1. Select "Scan Occlusion"
2. Once again, flip the scan tip over to its original position

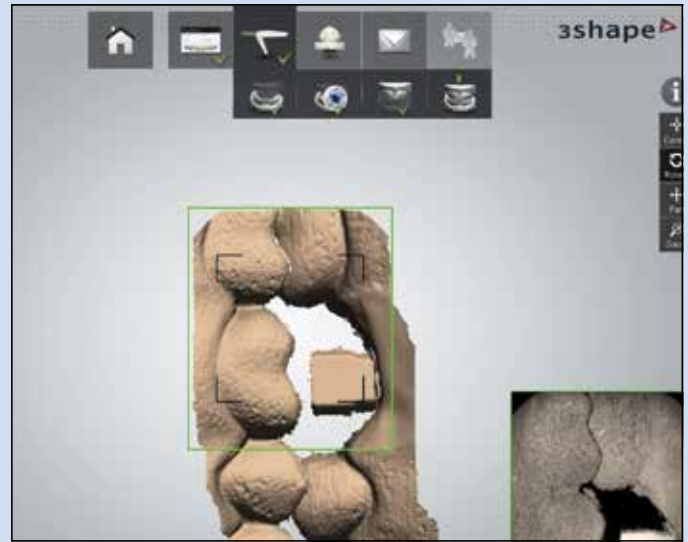


Scan the bite

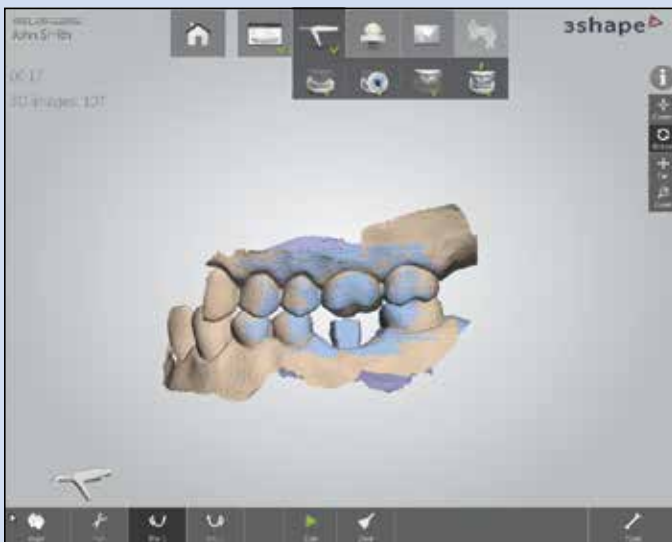


Bite has been established


How to Send Your 3Shape Files to Glidewell Laboratories

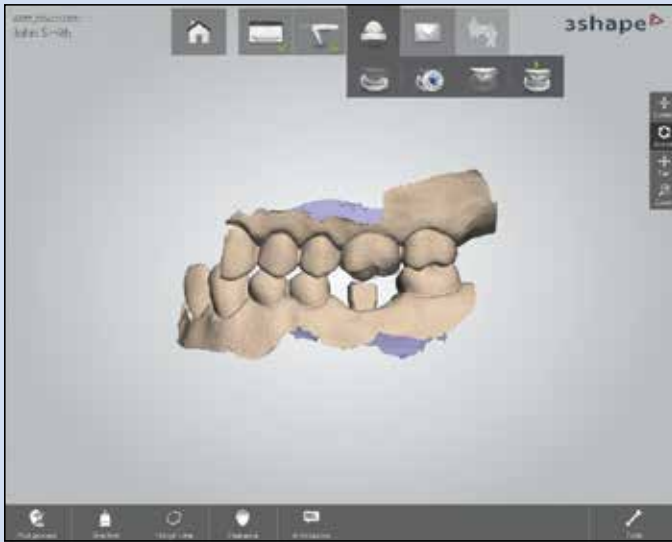

1. Press tooth icon
2. Press "Post Process"

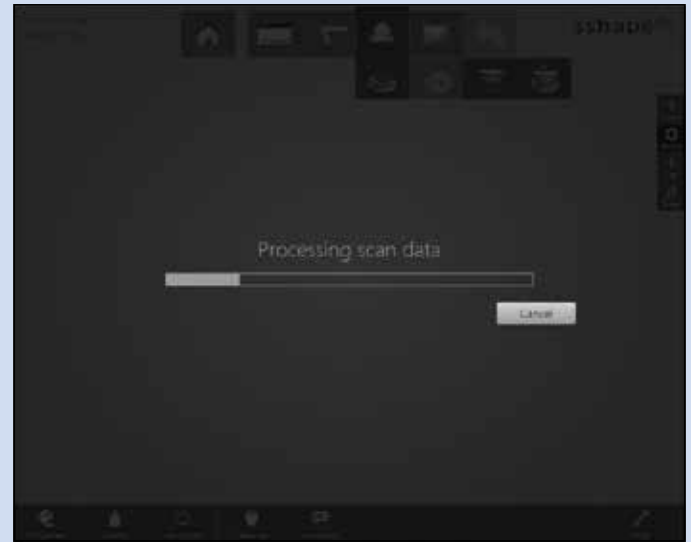

Processing scan data

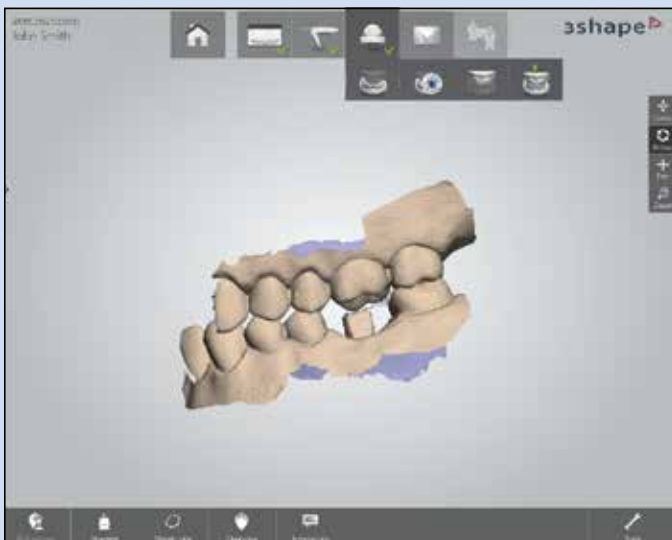

Press "Mail"

How to Send Your 3Shape Files to Glidewell Laboratories



Press "Send Order"



Order successfully sent



GLIDEWELL LABORATORIES

800-887-3580 • www.glidewelldental.com

GL-3388-0613