

BruxZir[®] *esthetic* NOW



PATIENT-PLEASING RESULTS



HIGH STRENGTH



NO OVEN NEEDED



< 40 MIN MILL TIME*

THE FULLY SINTERED, HIGHLY ESTHETIC ZIRCONIA FOR IN-OFFICE USE

BENEFITS

- Lifelike translucency engineered for the anterior
- Deliver restorations right from the mill with no post-processing required
- Preserve your patients' existing dentition with minimal preparation requirements

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184

glidewell.io[™]
IN-OFFICE SOLUTION



For this horizontal fracture on tooth #9, BruxZir® NOW Esthetic delivered the strength needed for a durable same-day restoration while blending seamlessly with the patient's natural smile.

SPECIFICATIONS

MATERIAL

- Pre-shaded, fully sintered zirconia
- Flexural strength: 870 MPa

IDEAL FOR


- Single-unit anterior restorations, including bicuspid
- Cases with limited interocclusal space
- Restorations with < 1 mm margin thickness

SHADES

- 16 VITA® Classical Shades: A1, A2, A3, A3.5, A4, B1, B2, B3, B4, C1, C2, C3, C4, D2, D3, and D4
- 2 Bleach Shades: BL1 and BL White

EASY PREPARATION

- A** 1.25 mm ideal reduction (0.7 mm minimum)
- B** Chamfer or modified shoulder margins preferred
- C** Axial walls must be convergent (avoid undercuts)
- D** Preparation should be cut in three planes
- E** To achieve optimal impression quality, gingival retraction is necessary for preparations with subgingival or equigingival margins

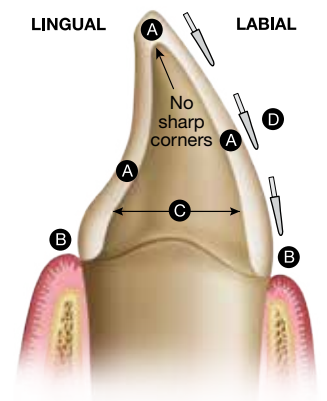

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184