

BruxZir[®] NOW

BRIDGE BLOCK



MAXIMUM STRENGTH



MINIMAL PREPARATION



NO OVEN NEEDED



< 2.5 HOURS MILL TIME*

THE ONLY FULLY SINTERED ZIRCONIA FOR IN-OFFICE USE

BENEFITS

- Deliver restorations right from the mill with no post-processing required
- Preserve your patients' existing dentition with minimal preparation requirements

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184

glidewell.io[™]
IN-OFFICE SOLUTION



This patient presented with failing restorations on #18–20, which were replaced with BruxZir® NOW Bridge Block restorations fabricated chairside with the glidewell.io™ In-Office Solution. The replacement BruxZir NOW restorations are esthetic, comfortable and strong enough to last the patient long into the future.

SPECIFICATIONS

MATERIAL

- Pre-shaded, fully sintered zirconia
- Flexural strength: > 1,000 MPa

IDEAL FOR

- Multiple-unit anterior and posterior restorations (average mill time: 2 hours, 30 minutes*)
- Cases with limited interocclusal space
- Restorations with < 1 mm margin thickness

SHADES


- 6 VITA® Classical Shades: A1, A2, A3, B1, C2 and D2

EASY PREPARATION


- 1.25 mm ideal reduction (0.7 mm minimum)
- Chamfer or modified shoulder margins preferred)
- Axial walls must be convergent (avoid undercuts)
- Preparation should be cut in three planes
- To achieve optimal impression quality, gingival retraction is necessary for preparations with subgingival or equigingival margins

- 1.0 mm ideal reduction (0.5 mm minimum)
- Chamfer or shoulder margins preferred. Feather-edge OK
- Axial walls must be convergent (avoid undercuts)
- Preparation should be cut in three planes
- To achieve optimal impression quality, gingival retraction is necessary for preparations with subgingival or equigingival margins)

ANTERIOR

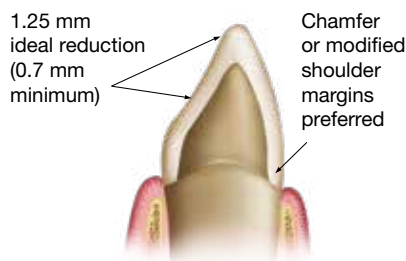

POSTERIOR

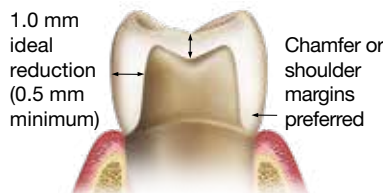

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184