


Obsidian[®] **NOW**



PATIENT-PLEASING
RESULTS



NO OVEN
NEEDED

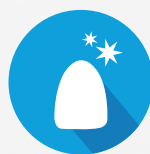

MINIMAL
PREPARATION



< 35 MIN
MILL TIME*

THE LITHIUM DISILICATE WITH EASY DELIVERY AND GREAT RESULTS

BENEFITS

- Deliver Obsidian[®] NOW right from the mill with no oven required
- Easy delivery — conventionally cemented or adhesively bonded
- Preserve your patients' existing dentition with minimal preparation requirements

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184

glidewell.io[™]
IN-OFFICE SOLUTION



Two anterior veneers on teeth #8 and #9 milled from Obsidian® NOW, delivering natural esthetics and a seamless blend with the surrounding dentition.

SPECIFICATIONS

MATERIAL

- Versatile monolithic glass-ceramic
- Flexural strength: 476 MPa

IDEAL FOR


- Single-unit anterior and posterior restorations
- Cases with limited interocclusal space
- Restorations with ≥ 1.2 mm thickness for circular and ≥ 1.5 mm for incisal/occlusal

SHADES

- Classical Shades: A1, A2, A3, B1

EASY PREPARATION

- Crown preparation should be tapered 6° along the long axis of the tooth
- Rounded internal with no sharp angles
- Margin thickness 1.0 mm minimum
- Palatal reduction 1.2 mm minimum
- Incisal edge reduction 1.5 mm minimum

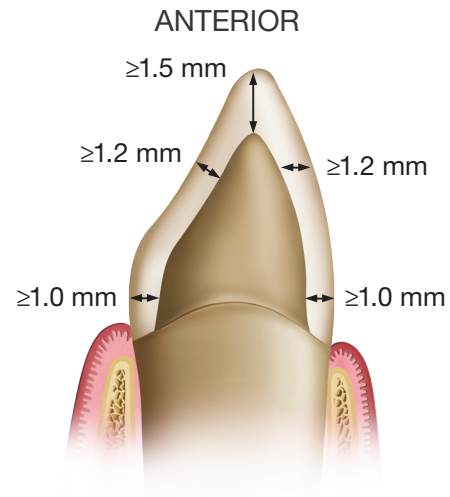

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184

*Milling times vary based on prescription and mill used. Average milling time is shown.

PK-5033462-061226