



Anterior

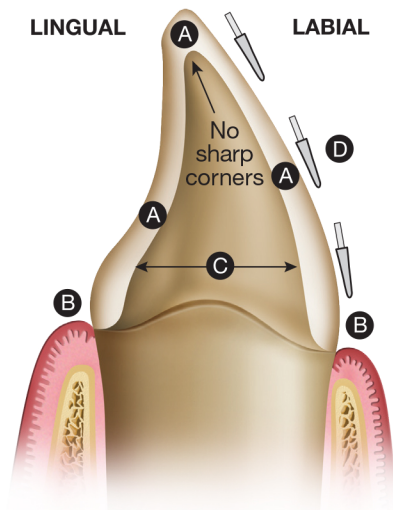

- A. 1.0 mm ideal reduction (0.5 mm minimum)
- B. Chamfer or modified shoulder margins preferred.
Feather-edge OK
- C. Axial walls must be convergent (avoid undercuts)
- D. Preparation should be cut in three planes
- E. To achieve optimal impression quality, gingival retraction is necessary for preparations with subgingival or equigingival margins

Posterior

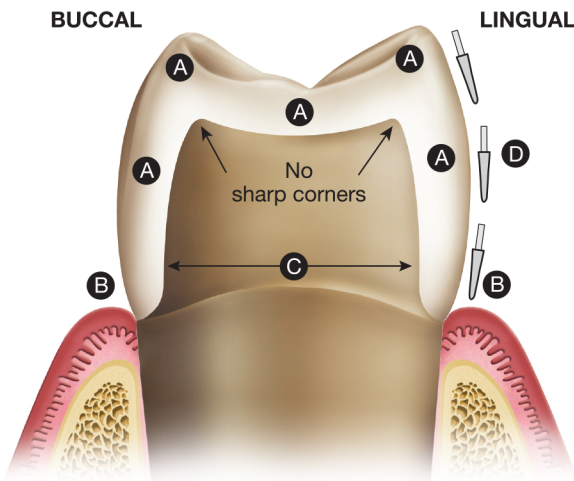

- A. 1.0 mm ideal reduction (0.5 mm minimum)
- B. Chamfer or modified shoulder margins preferred.
Feather-edge OK
- C. Axial walls must be convergent (avoid undercuts)
- D. Preparation should be cut in three planes
- E. To achieve optimal impression quality, gingival retraction is necessary for preparations with subgingival or equigingival margins