

CHOOSE THE SIMPLE, INTUITIVE, HANDS-OFF 3D PRINTING SOLUTION

fastprint.io™ 3D PRINTING SOLUTION




For more information about fastprint.io™
scan the QR code or visit:

glidewell.io/fastprint


844-949-7184




EASILY PRINT CRYSTAL-CLEAR APPLIANCES

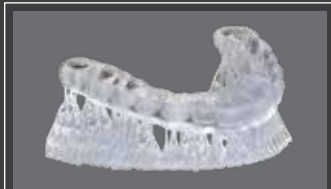

Appliances from other printers are cloudy and use excessive material

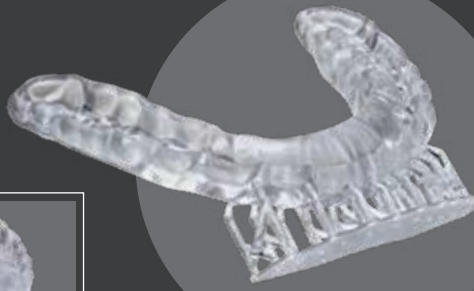

The fastprint.io™ uses LuxCreo Digital Polishing™ for crystal-clear appliances with less material waste

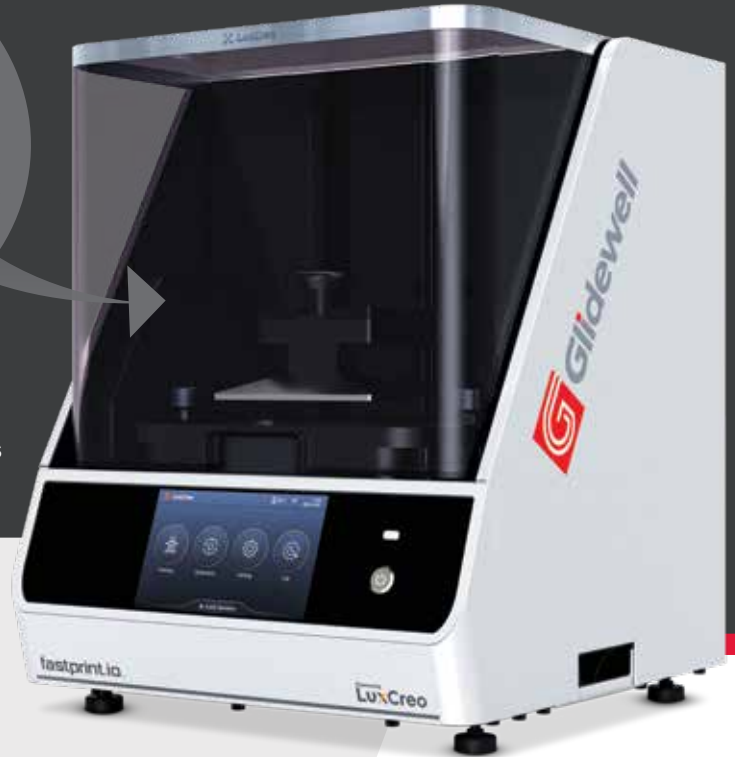

Use a range of Glidewell resins to easily print models and nightguards, plus a variety of LuxCreo-compatible materials to create crowns, retainers, and aligners.

- Use LuxCreo Digital Polishing™ to print clear appliances with no hand polishing needed
- Depend on AI-enhanced automation for a streamlined 3D printing workflow
- Elevate your clinical possibilities with a practical same-visit solution



COLLECT

The fastprint.io integrates seamlessly with the fastscan.io™ Scanning Solution and is compatible with all popular intraoral scanners on the market.



CONFIRM

Enjoy AI-driven print-ready designs and fast one-click staging, slicing, and printing with a seamless workflow that eliminates clunky software integrations.



CREATE

Print same-day aligners with the first FDA-class II, 510(k) cleared direct-print clear aligner workflow plus a range of other appliances for same-visit solutions.



COMPLETE

Minimal post-processing with LuxCreo Digital Polishing technology and an automated cleaning station allows for simple, same-visit delivery.

POWERED BY AI AND GLIDEWELL LAB SOFTWARE

The fastprint.io™ 3D Printing Solution utilizes the same AI-backed software used at Glidewell to create thousands of appliances every year, simplified for chairside use.

Confirm your AI-generated designs for nightguards or models in a single click, then send them to print. With this intuitive AI-backed program, there's no need to spend valuable time learning new software.

LuxCreo's LuxDental software continues the hands-off automated design process by automatically staging and slicing in preparation for print.




COMPLETE WITH AUTOMATED WASHING AND CURING STATIONS

LuxCreo's washing unit reduces the amount of labor in your workflow by enabling fully automated and unattended operations. With integrated heating, agitation, and ultrasonic functions, your printed parts come out with significantly less residual resin compared to traditional methods.

Your washing unit is also paired with a compact automated curing station designed for fast, reliable cure times. With a high intensity light that heats from above and below, there's no need to flip the devices during the curing stage.

Paired with LuxCreo's Digital Polishing™ technology which eliminates more than 85% of labor*, simplified 3D-printed same-visit solutions are now a reality.

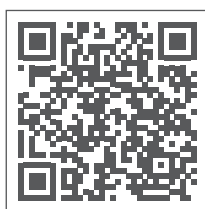


Scan the QR code to see fastprint.io in action.

*Based on data provided by LuxCreo.

TECHNICAL SPECIFICATIONS

Build Volume (XYZ)	6 x 3.4 x 7.9 in (153 x 86 x 200 mm)
Resolution	3840 x 2400 ppi
Printing Speed	Material dependent
3D Printer Dimensions (WxDxH)	18 x 20.8 x 22.9 in (457 x 528 x 582 mm)
3D Printer Dimensions with Door Open	18 x 20.8 x 28.7 in (457 x 528 x 729 mm)
3D Printer Weight	69.4 lbs (31.5 kg)
Wash Station Dimensions (WxDxH)	7.9 x 19.3 x 17 in (200 x 490 x 432 mm)
Wash Station Dimensions with Door Open	7.9 x 20.3 x 24.9 in (200 x 516 x 633 mm)
Wash Station Weight	38.6 lbs (17.5 kg)
Cure Station Dimensions (WxDxH)	8.1 x 11.5 x 11.4 in (205 x 292 x 290 mm)
Cure Station Dimensions with Door Open	8.1 x 18.1 x 11.4 in (205 x 460 x 290 mm)
Cure Station Weight	20.9 lbs (9.5 kg)
Warranty	12-month manufacturer warranty included Extended warranty options available

WORKSPACE REQUIREMENTS

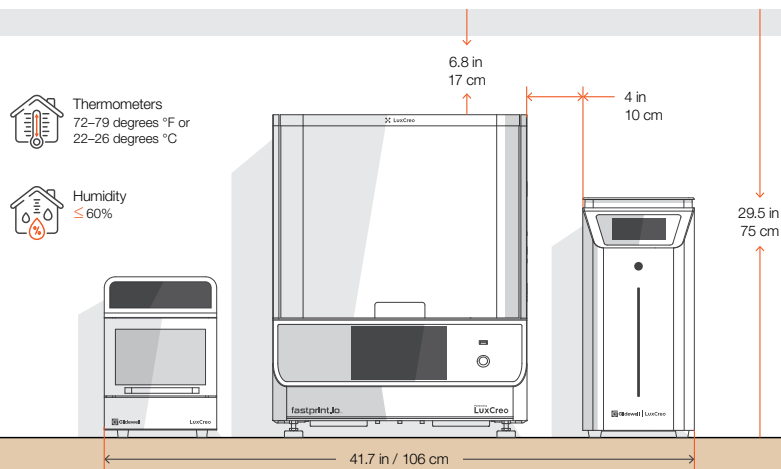

For more information regarding
workspace requirements
scan the QR code.

Ready to start printing your own clear appliances?
The choice is clear with the fastprint.io™ 3D Printing Solution.

Start printing with fastprint.io for only \$20,990*

Comes complete with washing and curing stations plus starter kit.
Contact us for our latest promotions and financing options!

*Price does not include shipping or applicable taxes.



For more information about fastprint.io
scan the QR code or visit:

glidewell.io/fastprint
844-949-7184

 **Glidewell**

 **LuxCreo**