

BruxZir[®] NOW



MAXIMUM STRENGTH



MINIMAL PREPARATION



NO OVEN NEEDED



< 40 MIN MILL TIME*

BRUXZIR NOW: THE FULLY SINTERED ZIRCONIA FOR IN-OFFICE USE

BENEFITS

- Deliver restorations right from the mill with no post-processing required
- Preserve your patients' existing dentition with minimal preparation requirements

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184

glidewell.ioTM
IN-OFFICE SOLUTION



This patient presented with failing restorations and recurrent caries on teeth #12-15. The replacement BruxZir® NOW restorations are esthetic, comfortable and strong enough to last the patient long into the future.

SPECIFICATIONS

MATERIAL

- Pre-shaded, fully sintered zirconia
- Flexural strength: > 1,000 MPa

IDEAL FOR


- Single-unit posterior restorations
- Cases with limited interocclusal space
- Restorations with $\geq 0.5\text{mm}$ for circular and $\geq 0.7\text{mm}$ for incisal/occlusal restorations

SHADES

- 16 VITA Classical Shades: (A1, A2, A3, A3.5, A4, B1, B2, B3, B4, C1, C2, C3, C4, D2, D3, D4)
- 3 Bleach Shades: Bleach White, Bleach 1, Bleach 3

EASY PREPARATION

- 0.7 mm incisal/occlusal reduction required (1.5 mm recommended)
- Chamfer or shoulder margins preferred; feather-edge OK
- All axial walls must be convergent

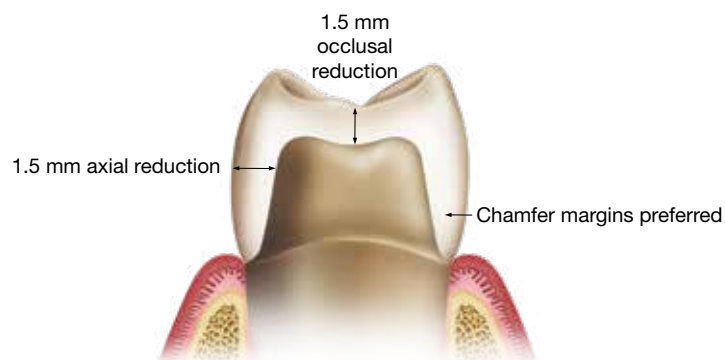

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184