

CAMouflage® *NOW*



RETAINS
SHADE



RESISTS
STAINS



NO OVEN
NEEDED



< 10 MIN
MILL TIME*

THE COMBINATION OF **SPEED, ESTHETICS AND DURABILITY**
YOU DESIRE

BENEFITS

- Strong survival performance among ceramic-polymer materials
- Fast delivery with no oven wait
- Easy to finish with high gloss retention
- Choose from a range of lifelike shades

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184

glidewell.io™
IN-OFFICE SOLUTION



A CAMouflage NOW overlay on #30 and onlay on #31 provided a minimally invasive single-visit solution.

SPECIFICATIONS

MATERIAL

- Approximately 80% ceramic-filled nanohybrid resin composite
- Flexural strength: 193 MPa

IDEAL FOR

- Single-unit anterior and posterior crowns
- Inlays and onlays
- Veneers


SHADES

- 7 VITA® Classical Shades, 5 in 2 translucencies: A1HT/LT, A2HT/LT, A3HT/LT, A3.5HT/LT, B1HT/LT, C2LT, D2LT
- 1 VITA Bleach Shade: BLLT


EASY PREPARATION

- Prepare cavity walls tapered 5–6° to the long axis.
- Wall thickness should be ≥ 1.5 mm.
- Isthmus thickness should be ≥ 1.5 mm.


ANTERIOR

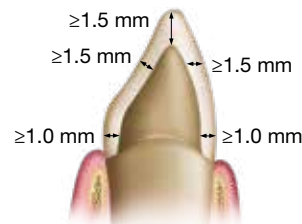


POSTERIOR

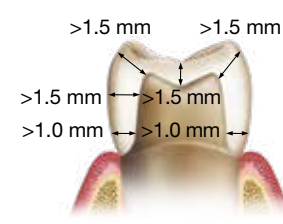


VENEERS

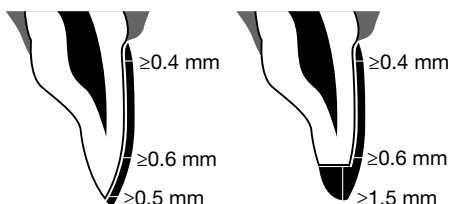

INLAY

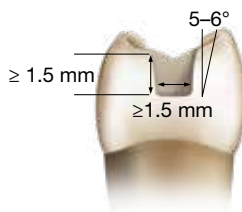

ONLAY

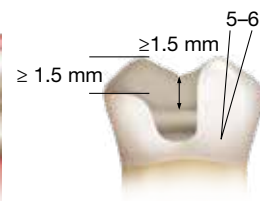

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184