



**SUPERIOR CHIP
RESISTANCE**



**PATIENT-PLEASING
RESULTS**



**WEAR-
FRIENDLY**



**< 20 MIN
MILL TIME***

THE LITHIUM SILICATE WITH EASY DELIVERY AND GREAT RESULTS

BENEFITS

- Provide patients with lab-quality Obsidian® crowns in just one appointment
- Select from bonding or cementing delivery methods to suit each case
- Choose from a wide range of lifelike shades that blend well with existing dentition

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184



This patient presented with failing PFMs and recurrent decay on teeth #7 and #8. Obsidian® crowns were delivered in just one visit. To improve overall symmetry, same-day Obsidian veneers were also placed on teeth #9 and #10.

SPECIFICATIONS

MATERIAL

- Lithium silicate ceramic
- Flexural strength: 385 MPa

IDEAL FOR


- Full-coverage crowns
- Inlays and onlays
- Veneers

SHADES

- 12 VITA® Classical Shades: A1, A2, A3, A3.5, B1, B2, B3, C1, C2, C3, D2, D3
- 2 Ivoclar™ Bleach Shades: BL1, BL4

EASY PREPARATION

- Tooth should be prepared with an expulsive 6° minimum angle taper along the long axis.
- Partial coverage restorations require wall thickness ≥ 1.5 mm
- Isthmus thickness ≥ 1.0 mm

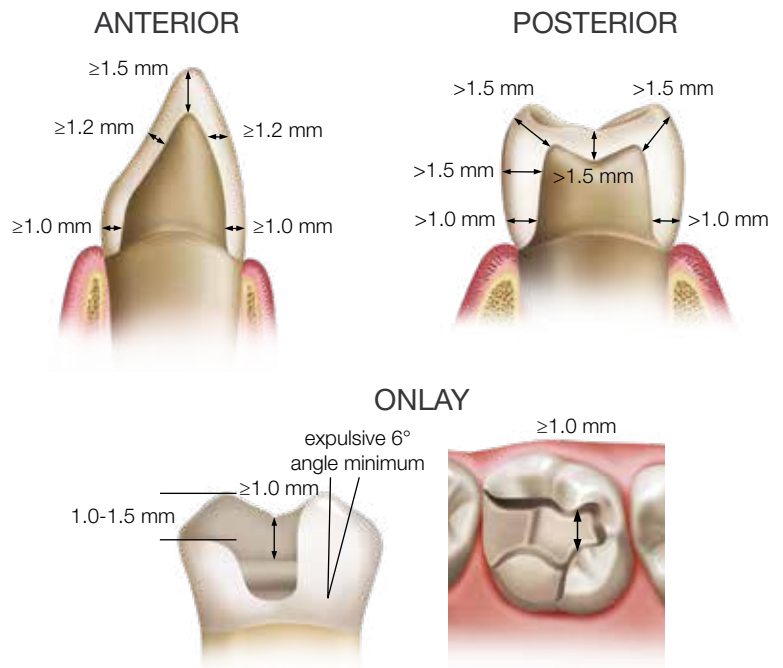

FOR MORE INFORMATION

glidewell.com/io-materials | 844-949-7184

*Milling times vary based on prescription and mill used. Average milling time is shown. VITA is a registered trademark of VITA Zahnfabrik. Ivoclar is a registered trademark of Ivoclar Vivadent.

PK-7294847-061226